


Course Instructor: Md. Nazmul Abedin Khan
Course Code: 109
Course Title: CHE

Student name: B M Shahria Alam
ID: 2021-3-60-016
Section: 5

Experiment-7: Determination of the acid dissociation constant (K_a) of a strong acid (hydrochloric acid, HCl).

Theory:

When an acid is dissolved in water, it breaks apart or dissociates to a slight extent. A proton from the acid is donated to a water molecule. The equations for the equilibrium and the equilibrium constant expression are as follows:

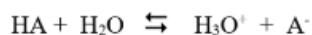

$$K_a = \frac{[H_3O^+][A^-]}{[HA]}$$

Where A^- represents the anion of the weak acid and the square brackets indicate molar concentrations of the species. The value of the equilibrium constant, K_a , indicates to what extent the reaction occurs. The greater the value of K_a , the stronger the acid, and the greater the amount of dissociation. HCl and Cl^- ions are conjugate acid-base pairs. A conjugate acid is a substance that has one more proton in its structure than its corresponding conjugate base. This combination also results from a mixture of a strong acid (HCl), and its salt (Cl^-). The equilibrium constant expression is:

$$K_a = [H_3O^+][Cl^-]/[HCl]$$

If a solution contains equal concentrations of HCl and Cl^- - these concentration terms cancel out in the above equation so that $K_a = [H_3O^+]$.

Procedure:

1. Fill the burette with 0.1M standardized NaOH solution.
2. Take exactly 50.0 mL of HCl (unknown concentration) in a beaker and record the pH of that acid solution.
3. Pour precisely 25.0 mL of the acid solution into an Erlenmeyer flask and add 2 drops (0.1 mL) of phenolphthalein solution.
4. Titrate this solution with NaOH solution by continuously adding small volumes (~0.2 mL or less) until slightly changing its color to a permanent pink. This process converts all the acid (HCl) in the flask into its conjugate base (Cl^-) according to the neutralization reaction ($OH^- + HA = H_2O + A^-$). Record the pH of the solution.
5. To calculate the K_a of the acid, follow any one of the following calculations (a or b).
 - (a) Pour the contents of the flask into the beaker and mix the solution. Measure the pH of this solution using both a pH meter. The pH is the pK_a of the acid. Calculate the value of K_a of the acid. (note: $pK_a = -\log K_a$).
 - (b) Calculate the final pH of the mixture of solution 1 and solution 2 mathematically to follow the equation:

$$C_1 V_1 + C_2 V_2 = C_f V_f$$

$$\text{Or, } C_f = (C_1 V_1 + C_2 V_2) / V_f$$

Where,

C1 : [H +] of solution 1; C 2 : [H +] of solution 2,

V1 : volume of solution 1 (25 mL of acid);

V2 : volume of solution 2 (25 mL of acid + added NaOH);

Vf : final volume of (solution 1 + solution 2)

Cf : final [H +] of (solution 1 + solution 2) =?


Therefore, the pH = -log [Cf]; The pH is the pKa of the acid.

Calculate the value of Ka of the acid. (note: pKa = -log K a).

Data:

Volume of NaOH, ml	pH of Acid
0	1
0.5	1.02
1	1.04
1.5	1.05
2	1.07
2.5	1.09
3	1.11
3.5	1.12
4	1.14
4.5	1.16
5	1.18
5.5	1.2
6	1.21
6.5	1.23
7	1.25
7.5	1.27
8	1.29
8.5	1.31
9	1.33
9.5	1.35
10	1.37
10.5	1.39
11	1.41
11.5	1.43
12	1.46
12.5	1.48
13	1.5
13.5	1.53
14	1.55
14.5	1.58
15	1.6
15.5	1.63
16	1.66
16.5	1.69
17	1.72

17.5	1.75
18	1.79
18.5	1.83
19	1.87
19.5	1.91
20	1.96
20.5	2.01
21	2.06
21.5	2.12
22	2.2
22.5	2.28
23	2.38
23.5	2.51
24	2.69
24.5	3
25	6.99
25.5	10.99

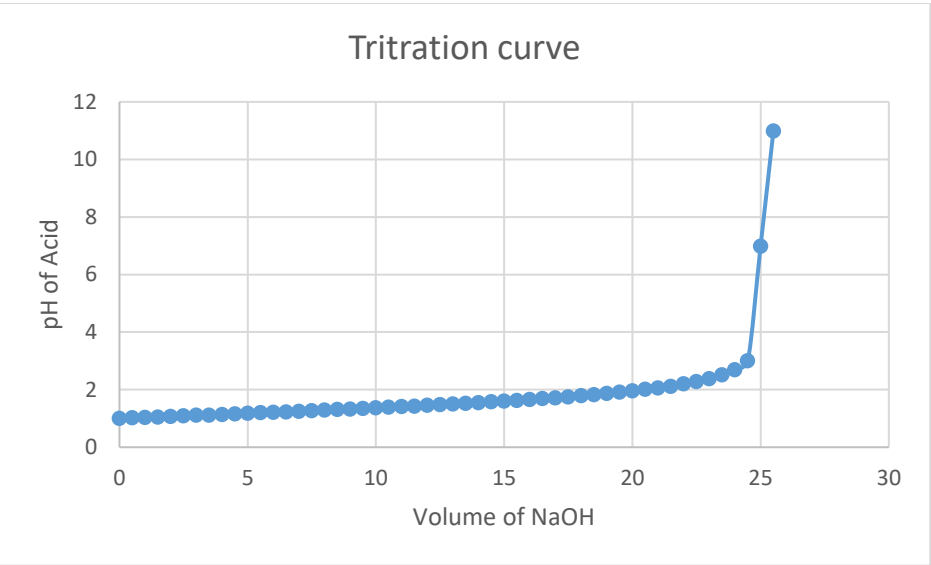

pH= pKa + log [salt]/ [acid]

pH of solution= 1.49
pH= pKa + log [25]/ [25]
pH= pKa + log [1]
pH= pKa + 0
pH= pKa= 1.49
Now, pKa= -logKa
Ka = 0.032359366

Virtual Lab File Edit View Help

EN Strong Acid Textbook Problems

Stockroom +

Information ≡

Name: 25ml (0.1 M HCl)
Volume: 0.0000 mL

Species (aq) Molarity

Display Absorbance


[Download Absorbance Table](#)

Temperature: 25.00°C


25.0 deg

pH: ∞


Workbench 1




0.1M NaOH
50.000 mL @ 25.0°C




25ml (0.1 M HCl)
0.0000 mL @ 25.0°C



Phenolphthalein
99.900 mL @ 25.0°C




Salt
75.600 mL @ 25.0°C



50 mL Burette
24.500 mL @ 25.0°C

1	Volume of NaOH, ml	pH of Acid
2	0	1
3	0.5	1.02
4	1	1.04
5	1.5	1.05
6	2	1.07
7	2.5	1.09
8	3	1.11
9	3.5	1.12
10	4	1.14
11	4.5	1.16
12	5	1.18
13	5.5	1.2
14	6	1.21
15	6.5	1.23
16	7	1.25
17	7.5	1.27
18	8	1.29
19	8.5	1.31
20	9	1.33
21	9.5	1.35
22	10	1.37
23	10.5	1.39
24	11	1.41
25	11.5	1.43
26	12	1.46
27	12.5	1.48
28	13	1.5
29	13.5	1.53
30	14	1.55
31	14.5	1.58
32	15	1.6
33	15.5	1.63
34	16	1.66
35	16.5	1.69
36	17	1.72
37	17.5	1.75
38	18	1.79
39	18.5	1.83
40	19	1.87
41	19.5	1.91
42	20	1.96
43	20.5	2.01
44	21	2.06
45	21.5	2.12
46	22	2.2
47	22.5	2.28
48	23	2.38
49	23.5	2.51
50	24	2.69
51	24.5	3
52	25	6.99
53	25.5	10.99

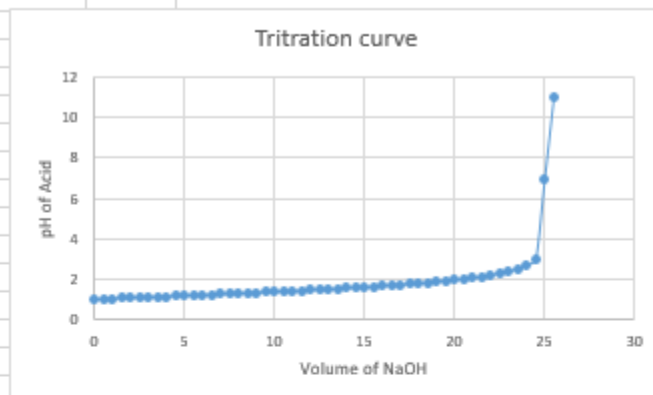

$$\text{pH} = \text{pKa} + \log \frac{[\text{salt}]}{[\text{acid}]}$$

$$\text{pH of solution} = 1.49$$

$$\text{pH} = \text{pKa} + \log \frac{[50]}{[50]}$$

$$\text{pH} = \text{pKa} + \log [1]$$

$$\text{pH} = \text{pKa} + 0$$

$$\text{pH} = \text{pKa} = 1.49$$

$$\text{Now, } \text{pKa} = -\log \text{Ka}$$

$$\text{Ka} = 0.032359366$$